




Alpha Medical Instruments LLC


BIPOLAR BALLOON CATHETER

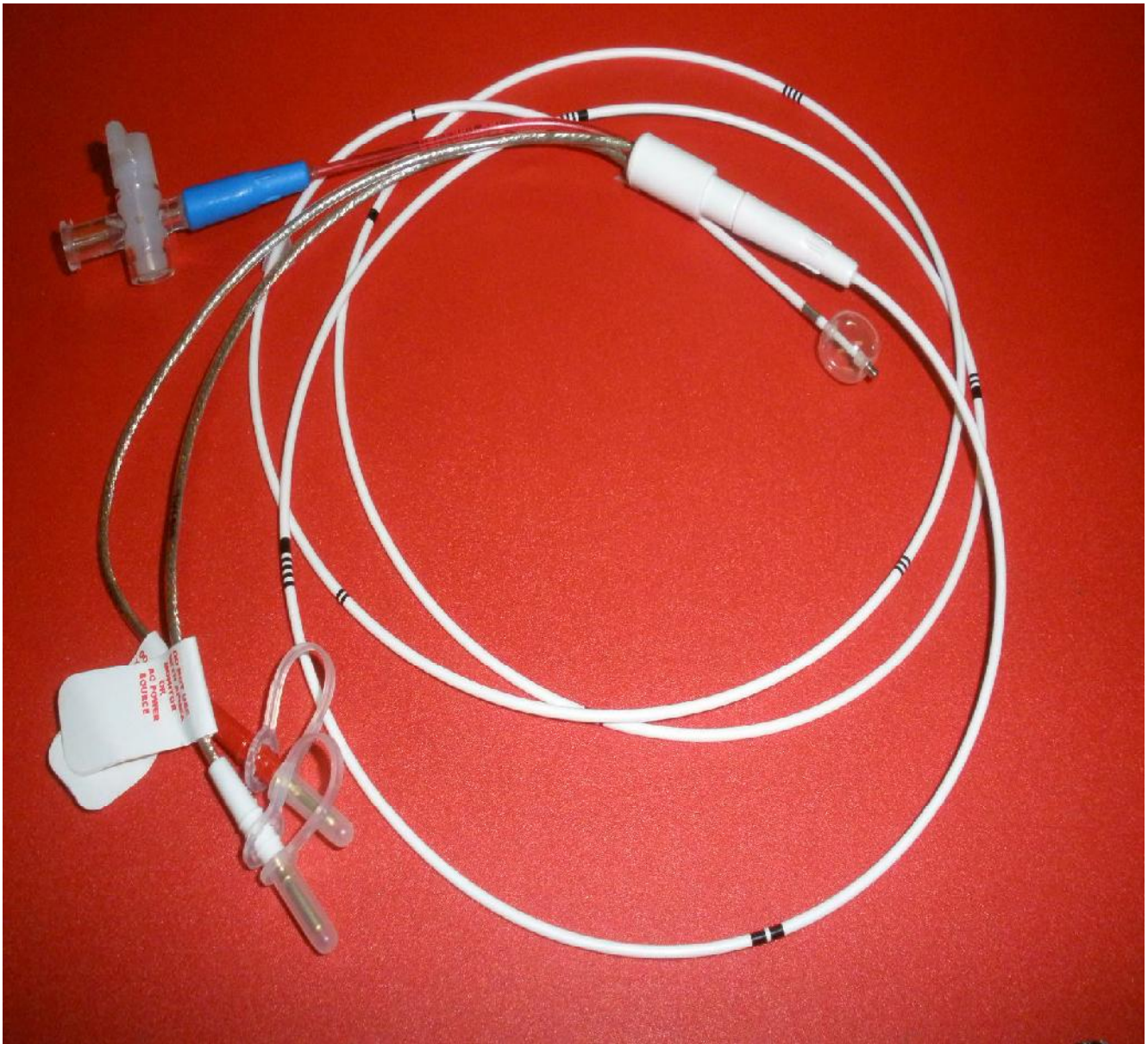

Bipolar Pacing Catheter

Series 600

Product Code:	Catheter Size:	Catheter Body Material:	Catheter Body Lumen:	Useable Length:	Inflated Balloon Diameter:	Measurement from Distal Tip Electrode Spacing:	Maximum Balloon Inflation:
652/1-90P	5 Fr.	PUR	2 L	90 cm	9 mm	1.0 cm	1.0 cc
652/1-110P	5 Fr.	PUR	2 L	110 cm	9 mm	1.0 cm	1.0 cc
652/2-90P	5 Fr.	PUR	2 L	90 cm	9 mm	2.5 cm	1.0 cc
652/2-110P	5 Fr.	PUR	2 L	110 cm	9 mm	2.5 cm	1.0 cc

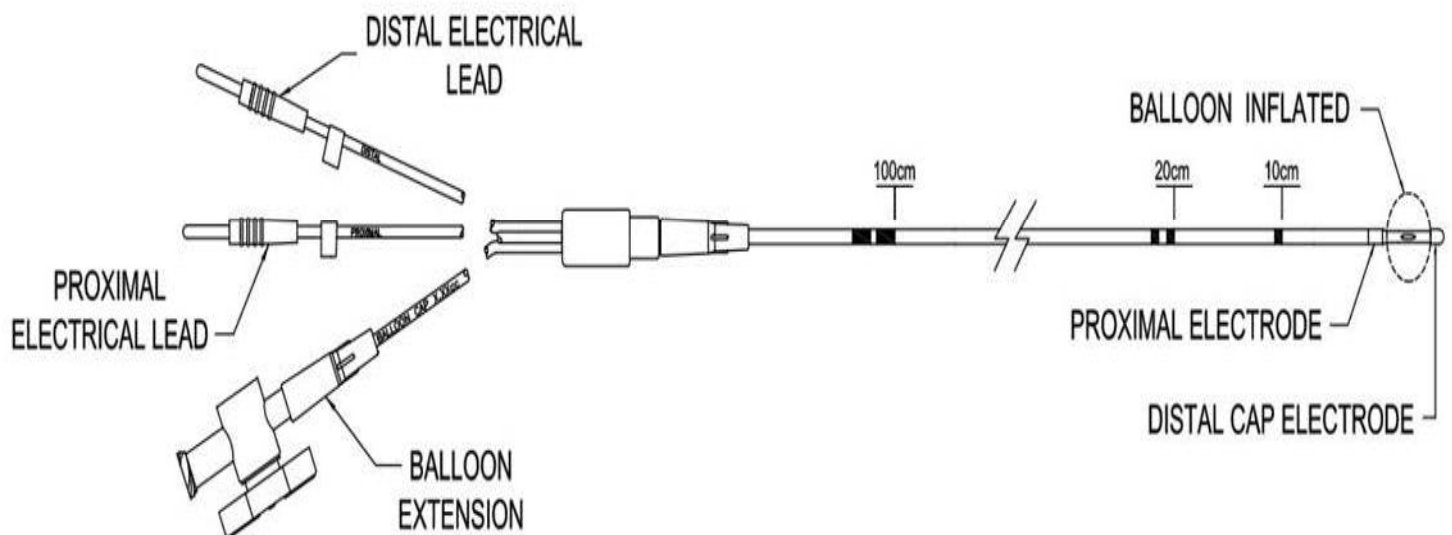

Description:

The Bipolar Balloon Pacing Catheter is made of extruded radiopaque polymeric tubing with two lumens and two stainless steel electrodes. Electrode spacing is 1.0 or 2.5 cm apart, depending upon which model you choose.

Packaging:

Each catheter is packaged in a sterile peel-pack tray

Catheters Length Marking:

Each catheter has clear, easy to identify increment markings every 10 cm along catheter body